



QUICK START GUIDE – How to Calibrate BRUNO MD6



READ BEFORE USE



1

NOTE: For Measuring Blood Glucose, calibration is Not Required. Calibration is ONLY Required for B-Ketone Strips.




3

Each box has a unique coding strip. Calibrate the meter with the code key from each new box. Discard the old key when using a new box of ketone strips.



5

Make sure the numbers on the key and box match.

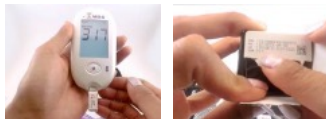

2

REQUIREMENTS CHECKLIST:

Before you do any testing, or handle sterile test strips, always ensure that you have washed your hands.

From the kit, you'll need:

A. Monitor, lancing device, and lancets.

B. Box of ketone strips (individually wrapped) & coding key (inside the box).



4

TO BEGIN CODING:

A. Insert the coding key.

B. If correct, you'll hear a "beep" & see CH briefly.

C. A number will appear on the display.



6

Congratulations! Your device is coded and ready!
Remove the code key and start using it.



VIDEO TUTORIALS



Please watch the how-to videos
before setting up your Bruno MD6



QUICK
START GUIDE



UNBOXING



LAB-ACCURATE
BLOOD GLUCOSE
& KETONE MONITORING



B-KETONES
TESTING &
SAMPLING



B-KETONES
CALIBRATION

Need help? Email customercare@brunopharma.com for assistance!